The First SUS Multiport Manifold on the Market

Nordson MEDICAL is proud to introduce the first all-in-one manifold for pharmaceutical and single-use biopharmaceutical fluid transfer applications. The revolutionary design of this patent-pending, compact fitting not only saves space, but also simplifies use and assembly. With five outlet ports available in one component, gone are the many leak points found in traditional stopcock and Y connector manifolds. Ideal for downstream processing, multiple bags can easily be filled from one large bag. Flow is consistent as it travels around the inner body of the part, maintaining a steady pressure. The large flow selector is effortless to turn, and the lock-in-place feature gives you assurance that you are utilizing the desired port.



Features & Benefits

- Can accommodate up to 5 fluid lines
- Flow restricted to one port at a time
- No pressure drop as fluid moves around & through the ports
- Intuitive push-to-turn flow selector
- Easier to use and fewer leak points than stopcock and Y connector manifolds
- Materials meet USP Class VI & ISO 10993 criteria
- Can be daisy chained for additional ports
- Arrows to provide a clear indication of how to turn the flow selector and which port is in use
- Etched data code for easy access to product resources and material information

Specifications


Materials:

Body & Flow Selector - Gamma Stable Polycarbonate Internal Spring - 316 Stainless Steel Gasket - Platinum-Cured Silicone

- Size: 3/8" (9.5 mm) tubing IDs
- Maximum Pressure: 15 psi (1.0 bar) liquid rated
- Barb Size:600 Series





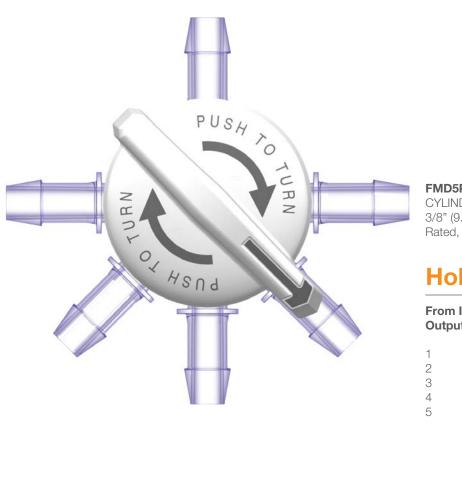


FMD5P670-001


CYLINDRAFlow[™] Manifold Director, 600 Series Barbs, 3/8" (9.5 mm) Flow, 1 Input, 5 Outputs, 15 psi Liquid Rated, -9024 Polycarbonate Material

Holdup Volume

From Inlet Port to	Holdup
Output Port Number	Volume
1	3.5 cc
2	3.9 cc
3	4.3 cc
4	4.7 cc
5	5.1 cc

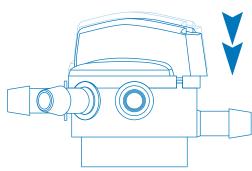

Push and Turn to Open Flow

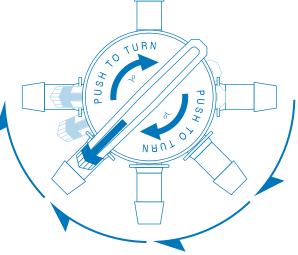

中国代理商:



捷阀流体技术(深圳)有限公司

- ◎ 深圳市龙岗区五和大道4012号 IOT物联网产业园一栋
- \$ +86 755 8452 8460/80
- enquiry@gefa-ultrapure.com.cn
- www.gefa-ultrapure.com.cn



